



We are a French company based in Bordeaux since 2014.
Our team specializes in the manufacture and distribution of innovative mobility aids that are design, lightweight, practical and easy to transport.

Our mission is to transform the lives of people with reduced mobility by offering them products that **combine comfort**, **safety and independence**.

Some of our products:



Manufacturer of Yomper+, power add-on for manual wheelchair users.



Protects arms and shoulders of the users by give them a smooth, intuitive and easy to handle motorization.





French distributor of folding electric wheelchairs SmartChair.



Regain your independence with a range of folding electric wheelchairs that are easy to transport and can be adjusted to the user's needs.



# **ENJC**

Manufacturer of Enjo, power add-on for manual wheelchair caregiver.



Relieve caregivers with a practical, intuitive and lightweight pushing aid.





Up to 35 KM

Only 5.8 KG Max 150 KG

→ Rigid and foldable wheelchairs

★ Attendant control option

→ Slopes up to 15%





#### **EASY TO DRIVE**

#### Intuitive and reassuring

The combination of its innovative positioning in-line with the rear axle of your wheelchair and its freewheeling capability allows you to perfectly maintain the dynamics of your wheelchair.

- + I steer my wheelchair with the handrims.
- + I brake my wheelchair with the handrims (indoor mode).
- + I pilot thanks to a simple wireless controller.
- \*I start smoothly and safely.

Start

Change the speed

Power off

2 taps to start

With the selector speed 4

Inside mode Grab handrims
Outside mode 1 tap

### **EVERYWHERE WITH YOU**

#### Small but powerful

Thanks to a powerful 250W motor capable of speeds up to 10km/h, it facilitates your daily trips. Spare your shoulders and arms even when traveling in difficult conditions (slopes, grass, gravel...).

- + Compact and lightweight (5.8kg), it will be your travel companion (carrying case included).
- + Energy efficient, it has a range of up to 35km.
- + Two small wheels (at the front of your Yomper) and raised handle allow you to easily attach and/or detach it from your wheelchair.

International Air Transport Association (IATA) approved

# ips. vel...).



Keep your wheelchair

The Yomper+ is compatible with most of rigid or foldable wheelchairs. The folding chair remains foldable thanks to its retractable support!

- ► I can configure the Yomper+ according to my strength and my mobility requirements.
- ►The 5 speeds per mode (indoor or outdoor) allows me to adapt my speed to my environment.





You can connect a third-person control to the Yomper+ with a cable. The speed is controlled with a simple push button without any loss of performance.

An attendant can drive and control the wheelchair while walking behind it, without any effort.

As the motor is positioned in the middle of the wheelchair, there is no offset to the rear. Your caregiver will not have any discomfort at the level of his feet to come and push.







Up to 23 KM

Only 4 KG Max 150 KG

→ Fit under all manual wheelchairs → Easy to install and remove → Slopes up to 15%







#### RELIEVE YOUR CAREGIVER

#### Intuitive and reassuring

Pushing a wheelchair can be hard, especially when faced with slopes or long distances.

You can guickly experience fatigue, muscle pain and stress, which can limit you.

With our Enjo motorization, getting out is much easier.


#### FIT UNDER ALL MANUAL WHEELCHAIRS

#### Designed to be compact

We can fit manual fixed frame wheelchair, foldable ones and even comfort wheelchair which are really heavy to push!

Enjo is designed for wheelchairs with drumbrakes, or for users that don't need breaks.



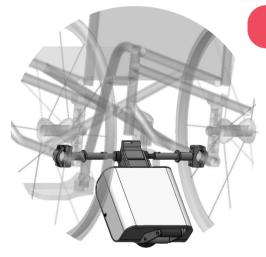

## EASY TO INSTALL AND REMOVE

#### Light and adjustable

You can install the Enjo in the blink of an eye.

Only 4kg to move, and you simply need to clip the motor to the attachment mounted on your wheelchair. And off you go!

The Enjo is design to be easily transfer from one chair to another in a matter of seconds.

This is a major advantage, particularly for rehabilitation centers, hospitals, retirement homes and cultural venues.

## CENTRAL POSITION IN THE WHEELCHAIR AXIS

#### Easy to handle

Benefiting from our 10 years experience with the Yomper, Enjo inherited the central axle position.

You keep a total maneuverability of the wheelchair even with the motor switched off.

This position, with no offset at the rear of the wheelchair, means that the user can push the wheelchair without any obstruction to his or her feet.



# They talk about it better than we do



I've been in a wheelchair since birth, and I've always been very sporty and independent. The Yomper seemed obvious to me because it's very discreet, light and practical, and it means I don't have to strain my arms as much, but I'm completely autonomous in everything. I can take it everywhere, I can use it everywhere.

It allows me to do a lot of things today, more than I did before.

#### ÉLODIE



I had an MVA and have been a paraplegic for 25 years.

I turned to the Yomper for its ease of storage, so I could move around by myself, install it on my own under my wheelchair and put it in the trunk of my car.

The Yomper has also enabled me to do more activities with my 12-year-old daughter.

It's really helped me to be able to do things I might not have been able to do before, or with more difficulty.

# SÉBASTIEN





I discovered Enjo at a trade show and was delighted by the idea of having a motorization that could be adapted to all the wheelchairs in the rehabilitation center where I work. This innovation will bring relief to hospital staff, elderly couples and carers in public places like airports, museums and so on.

The future looks bright for this little motor, which is also superbly designed.



